

DVMAX PACS Link

Digital Radiography Integrated Solution



A Powerful DICOM Solution

DVMAX PACS link is a newly integrated digital radiography feature in DVMAX Practice. DVMAX users can now **request digital radiography**, and then **view resulting medical images** (studies) right from the patient record within the DVMAX software. In addition, users can **launch a selected study** directly from DVMAX into the **DICOM viewer of their choice** for detailed analysis and manipulation.

DVMAX and SecureRad

DVMAX PACS link is the result of a joint effort between DVMAX and SecureRad, a PACS solution provider. This partnership allows DVMAX to offer a PACS system that is much less expensive and more flexible than other systems on the market. See www.securerad.com for more information on the PACS system.

Request digital radiography in DVMAX – integrated MWL (modality work list)	✓
View resulting medical images in DVMAX – within the medical history	✓
Open medical images into the DICOM viewer of your choice from the Patient record	✓
Distributed viewing on multiple workstations - no additional licensing fees	✓
Offsite (“cloud”) storage included – provides easy offsite access to your records	✓
Fast transmission speeds - download speed of several seconds	✓
Flexible storage solutions - power and growth potential	✓
Less expensive than other systems on the market - Lower cost per study	✓

Call 800 877 9221 x 1 today!


Medical Image Viewing

DVMAX and the SecureRad PACS software work together to provide full DICOM medical image viewing. Because the imaging software is based on a progressive open source product, and DICOM industry standards, it is constantly being enhanced to deliver the fastest transmission speeds possible. The system also includes remote “cloud” storage and access to your studies at no additional cost.

The imaging software also provides an integrated MWL (modality work list) that takes image requests from DVMAX and passes them on to your modalities (the medical imaging equipment) such as digital X-ray, Ultrasound, MRI and CT. As well, this software allows for distributed viewing on any DVMAX workstation with no additional licensing fees.

How it Works

DVMAX sends patient data to the PACS Link that turns the data into the proper modality work list (MWL) data format. The PACS then take the MWL Patient data and sends it to the appropriate Modality for processing. The resulting images/Studies are sent to the PACS for storage. The Images/Studes can then be retrieved and viewed in DVMAX or your choice of DICOM Viewer.

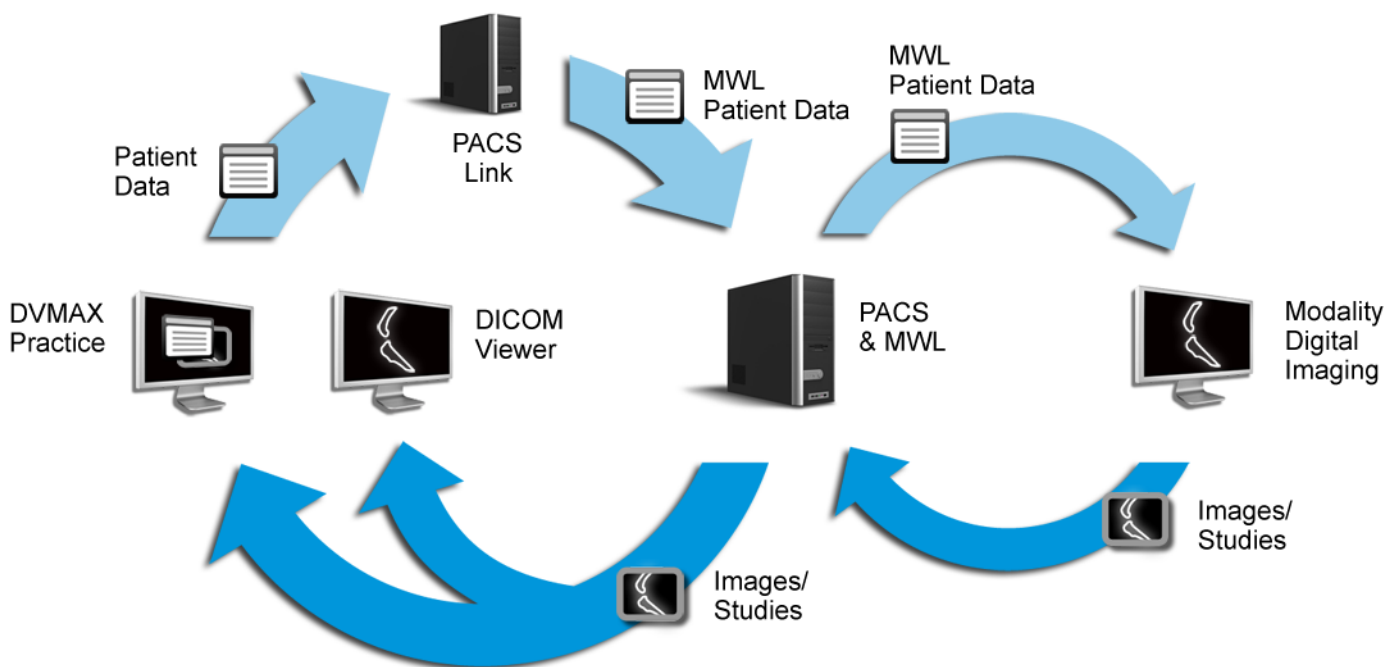

To see how a DVMAX user is taking advantage of the SecureRad PACS go to the following link:
<http://www.apple.com/science/profiles/avmi/index2.html>

Call 800 877 9221 x 1 today!